AL7010

AL7010 is designed to measure alcohol levels within the human breath.

The AL7010 displays clear blood alcohol concentration (BAC) results on the white LCD backlight screen. The accuracy of this device has been established at a Blood Alcohol Concentration of 0.000 to 0.400%BAC.



SPECIFICATIONS

1. Indication of B.A.C.: 0.000 to 0.400%BAC or 0.000 to 4.000% or 0.000 to 2.000mg/ ℓ BrAC by LCD display

2. Sensor: Oxide semiconductor type alcohol sensor

3. Accuracy: ± 0.01 %BAC up to 0.05%BAC

4. Mouthpiece: use

5. Power: Two 1.5V AA alkaline batteries (included)


6. Weights: 98g (Including Batteries)

7. Dimension: 105 (height) X 50 (width) X 18 (thickness) mm

8. Operation: 0 to 40 $^{\circ}$ C 9. Available color: Black

Contents

1x Breathalyzer, 5x Mouthpiece, 2x AA size Alkaline Batteries, 1x Zipper case, 1x Manual, 1x Hand strap

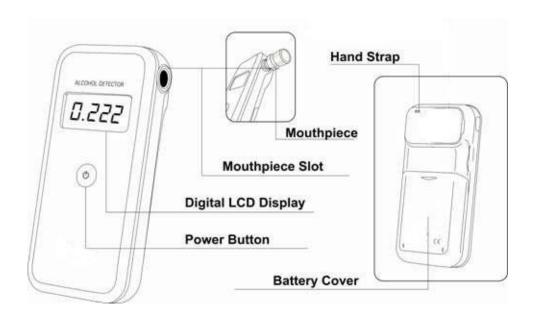

Sentech Korea Corp. Breath Alcohol Tester





Semiconductor sensor

Mouthpiece

Warning Display



If you see "Flow" in the display, the breath sample was not strong enough for measurement. Please repower the device and make sure the breath is strong and steady during sampling.

If you test over 300 tests, the AL7010 will display "CAL" sign for 15 seconds and then warm up for testing. The "CAL" sign will show when calibration is required.

If you see "bat" in the display, the battery power is low. Please replace the batteries.

Please contact

21-6, Jimokro 75-gil, Paju-Si, Gyunggi-Do, Korea 10880

TEL: 82-31-8071-4400 FAX: 82-31-8071-4411

URL: www.sentechkorea.com

E-mail: sentech@sentechkorea.com